




## PRODUCT INSTRUCTION SHEET

**PRODUCT CATEGORY:** Sensonics' Taste Strip Tests

**PRODUCT NAME:** Self-Administered Waterless Empirical Taste Test (WETT®-SA53)

### **DESCRIPTION:**

The Brief Self-Administered Waterless Empirical Taste Test™ (WETT®-SA53) is an easy and inexpensive method for quantitatively evaluating taste function. Delivery of tastants is via disposable plastic strips with proprietary pads of monomer cellulose that are embedded with either sucrose (sweet), citric acid (sour), sodium chloride (salty), caffeine (bitter), monosodium glutamate (umami), or blank tastants. No liquid rinsing is required. The test strips are located in four pouches or packs from which they can be individually removed for testing. The first pack contains strips numbered from 1 to 14, the second 15 to 27, the third 29 to 42, and the fourth 43 to 53.

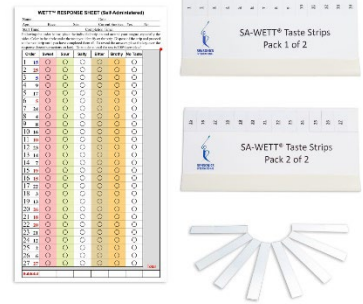

### **INSTRUCTIONS:**

- The WETT® should ideally be administered in a quiet test environment with no distractions.
- The WETT-SA53 should be placed on a table or appropriate workspace. After it's use, each taste strip should be discarded into a trash receptacle such as a wastebasket.
- The subject (or patient) should wear unscented latex gloves during testing. The taste strips should not be touched without gloves so as not to contaminate the tastants.
- Place each strip, one at a time, beginning with strip number 1 on and around your tongue, especially the sides, for 10 to 15 seconds. Using the provided response sheet, in the circle under the taste you identify on the strip as sweet, sour, salty, bitter, brothy (soupy, savory taste), or no taste at all.
- Dispose of the strip and proceed to the next strip until you have completed all strips.
- To reveal the answers, place the answer key over the response sheet.
- If the examiner does not want the subject to self-score the test, then the answer key should not be provided with the kit and the subject should not be given this element of the instruction.
- For most normal subjects, the administration test time is 5-10 minutes.



30°C  
Upper limit storage  
temperature



Keep dry



**Sensonics International**  
411 South Black Horse Pike  
Haddon Heights, NJ 08035  
United States  
Contact no. (856)-547-7702

MANUFACTURER [www.sensonics.com](http://www.sensonics.com)



**MedEnvoy Global B.V.**

Prinses Margrietplantsoen 33 - Suite 123  
2595 AM The Hague  
The Netherlands

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the member state